

NasoPore® FD

Fast Degrading Bioresorbable Nasal Dressing

2x faster degradation

With NasoPore FD nasal dressing, you'll get the same trusted benefits of previous NasoPore dressings, now in a fast degrading formula. It rapidly clears in an average of 5-7 days which is two times faster than NasoPore Firm dressing.




Support you need clearance you want

Your market-leading nasal dressing just got even better. Now with an average fragmentation time of 5-7 days, NasoPore FD nasal dressing widens your therapeutic options so you'll have a nasal dressing that's precisely right for each type of patient or procedure.

- 100% synthetic proprietary polymer
- 5-7 day average fragmentation time
- 36-48 hours tissue support/shape integrity
- Absorbs an average of 20x its weight (+/-8) without swelling in size
- Maintains moisture and clears through natural mucosal flow*
- No mixing, moistening, injection syringe or other prep needed


* Degradation may vary based on patient and environmental factors, amount of material inserted, and the type and manner of nasal rinsing

Fragmentation Time¹

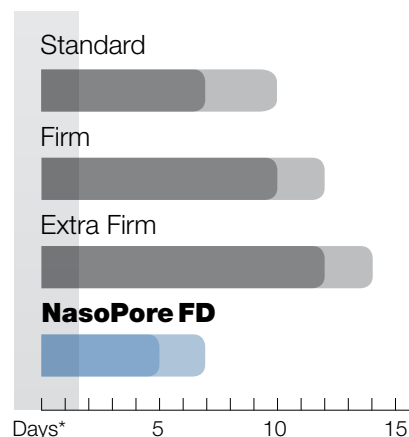

- ▲ Fragmenting in an average of 5-7 days, the dressing typically drains prior to the first post-op appointment. This provides a clear view to evaluate progress and eliminates manual removal which has potential for tissue injury or pain.

Compression Strength/Stenting¹

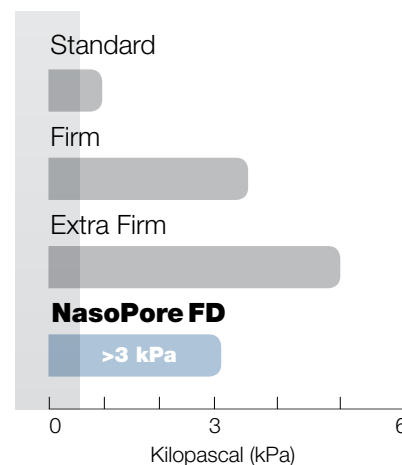

- ▲ The dressing's firm compression strength separates and stabilizes vital tissue structures such as the middle turbinate for 36-48 hours, thus preventing adhesions and minimizing edema.

Ordering is easy

Part Number: 5400-020-108 (8 per pack; 8cm size)

All variations of our bioresorbable solutions are conveniently available through a single source. To order, call Customer Service at **800-253-3210** or contact your Stryker NSE sales representative. For more details visit www.stryker.com/bioresorbables



References

1. *In vitro* data

The information presented in this brochure is intended to demonstrate the breadth of Stryker product offerings. Always refer to the package insert, product label and/or user instructions before using any Stryker product. Products may not be available in all markets. Product availability is subject to the regulatory or medical practices that govern individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Products referenced with TM designation are trademarks of Stryker. Products referenced with ® designation are registered trademarks of Stryker.

 Literature Number: 9100-003-407 Rev None G55

See NasoPore FD Instructions for Use (IFU) for complete details.

Copyright © 2015 Stryker
Printed in USA

Stryker Neuro Spine ENT
4100 East Milham Avenue
Kalamazoo, MI 49001 USA
t: 269 323 7700 f: 800 999 3811
toll free: 800 253 3210

www.stryker.com/NSE